

THE ASSORTMENT

BLADES & SCALPELS

The choice is yours!

Individual requirements, areas of application and personal habits play as much a role in the choice of the right scalpel or blade as do the product's technical specifications. These various factors form the foundation of our assortment of blades and scalpels.

Disposable scalpels

All HSO scalpel blades are made of **Swedish stainless steel**. The plastic handles are **ergonomically shaped** and come with a **metric scale**. The disposable safety scalpels have an additional retractable sleeve.

Reusable scalpel handles

In response to Directive 2010/32/EU on the prevention of sharp injuries in the hospital and healthcare sector, we also offer **ZEUS safety scalpel handles N°3 and N°4** to complement our safety scalpels. The handles, which are **reusable after proper sterilisation**, are made of stainless steel with an additional Teflon coating. The **uncoated metal protective slider** is pushed over the blade, like the mechanism of a utility knife, where it clicks into place to cover the blade completely and minimise the risk of injury. The blade is just as easy to unlock. After use, a lever in the metal protective slider helps to prevent accidental contact with the blade during its disposal.

Blades

The HSO range of blades fit N°3 and N°4 handles. The **standard blades** for scalpel handle N°3 (blades nos. 10–15SM) and for scalpel handle N°4 (blades nos. 18–36) are made of carbon steel for disposable use.

Herkules blades for scalpel handle N°3 (blades nos. 15SM, 15, 11 and 10) and handle N°4 (blades nos. 20, 21) are also made of carbon steel but their reinforced spines make them stronger and more break-resistant than conventional blades.


www.hso-solutions.at



SCALPELS BLADES



If you have any questions, please
contact your personal product specialist:

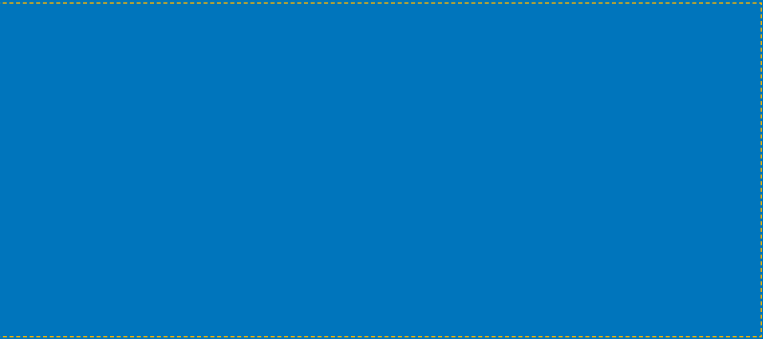

SALES IN AUSTRIA
Rudolf Heintel Gesellschaft m.b.H.
Erdbergstraße 166
1030 Wien, Österreich


T +43 (1) 403 89 56-0
F +43 (1) 403 89 56-50
medizintechnik@heintel.at
www.heintel.at

Photo:
© IStock Photo extravagantni
© IStock Photo -VICTOR-
© IStock Photo Cybernesco

Subject to errors and omissions.
Version 1.1 | Last updated: May 2021



ANY WAY YOU SLICE IT

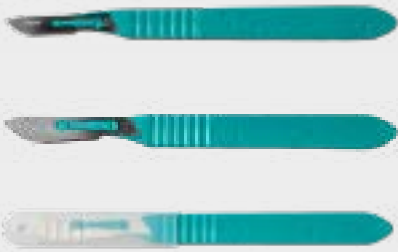

Blade no.	Order no.
10	HEIS10
11	HEIS11
11P	HEIS11P
12	HEIS12
12D	HEIS12D
15	HEIS15
15C	HEIS15C
15SM	HEIS15SM
20	HEIS20
21	HEIS21
22	HEIS22
23	HEIS23
24	HEIS24
25	HEIS25



Blade no.	Order no.
10	HSOSS10
11	HSOSS11
15	HSOSS11P
20	HSOSS20
21	HSOSS21
22	HSOSS22
23	HSOSS23
24	HSOSS24
25	HSOSS25

SAFETY-STITCH CUTTER



Blade no.	Order no.
Stitch Cutter	HSOSS-STITCH


Before using safety scalpels for the first time, it may be necessary to receive training by a certified product specialist. Please consult your national medical devices legislation.



Size	Order no.
3	HEISSG3
3L	HEISSG3L
4	HEISSG4
4L	HEISSG4L

Before using safety scalpels for the first time, it may be necessary to receive training by a certified product specialist. Please consult your national medical devices legislation.


N°3

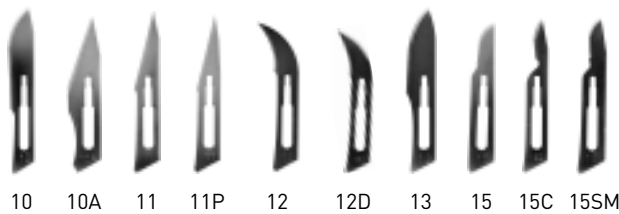

Blade no.	Order no.
10	HEIK10
10A	HEIK10A
11	HEIK11
11P	HEIK11P
12	HEIK12

Blade no.	Order no.
12D	HEIK12D
13	HEIK13
15	HEIK15
15C	HEIK15C
15SM	HEIK15SM


N°4

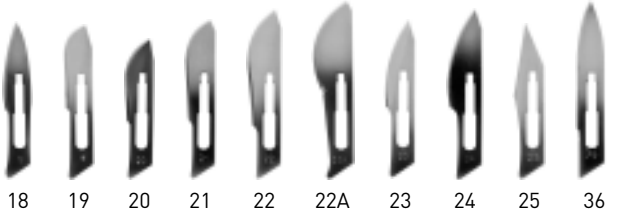

Blade no.	Order no.
18	HEIK18
19	HEIK19
20	HEIK20
21	HEIK21
22	HEIK22

Blade no.	Order no.
22A	HEIK22A
23	HEIK23
24	HEIK24
25	HEIK25
36	HEIK36

N°3

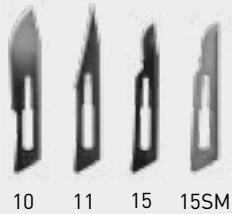

Blade no.	Order no.
10	HEIKH10
11	HEIKH11
15	HEIKH15
15SM	HEIKH15SM

N°4



Blade no.	Order no.
20	HEIKH20
21	HEIKH21
22	HEIKH22



Blade no.	Order no.
stitch cutter standard	HSOSTCS
stitch cutter midi	HSOSTCM
stitch cutter long	HSOSTCL